

## Chest Drainage Systems

Code	Description
NC-5512	2300ml System, High Capacity o
NC-3311	1150ml System, Compact

Qty/Box: 5

Set Contents

- 2300/1150ml Chest Drainage System
- 2/1 Tubing Set with Antikink Springs
- 2 Bed Hangers
- Filling Funnel

## Single Bottle System

Code	Description
NE-1101	Thorametrix Single Chest Drainage System

Qty/Box: 20

Set Contents

- 1100ml Single Bottle System
- 1 Tubing Set with Antikink Springs
- 1 Bed Hanger

## Silicone and PVC Radiopaque Thoracic Catheters

PVC		
Code	Description	Fr. Size
NM-1160	Straight	16
NM-1200	Straight	20
NM-1240	Straight	24
NM-1280	Straight	28
NM-1320	Straight	32
NM-1360	Straight	36
NM-1400	Straight	40
NM-2160	Right Angle	16
NM-2200	Right Angle	20
NM-2240	Right Angle	24
NM-2280	Right Angle	28
NM-2320	Right Angle	32
NM-2360	Right Angle	36
NM-2400	Right Angle	40

Qty/Box: 50

Silicone		
Code	Fr. Size	Length
NL-1620	16	70cm
NL-1820	18	70cm
NL-2020	20	70cm
NL-2420	24	70cm
NL-2820	28	70cm
NL-3220	32	70cm
NL-3620	36	70cm
MW-3190	19Fr Round	110cm
MW-3240	24Fr Round	110cm

Qty/Box: 50, MW-XXX0: 25

## Trocar Catheters

Code	Fr. Size	Catheter Length	Code	Fr. Size	Catheter Length
NN-0008	8	20cm	NN-0020	20	39cm
NN-0010	10	23cm	NN-0024	24	39cm
NN-0012	12	23cm	NN-0028	28	39cm
NN-0016	16	23cm	NN-0032	32	39cm

Qty/Box: 20


## Pneumothorax Sets

Code	Description
NP-1000	Non Seldinger Type Pneumothorax Set
NP-1005	Seldinger Type Pneumothorax Set
NP-1001	Seldinger Type Pneumothorax Set with Drainage Bag

Qty/Box: 10

# THORAMETRIX™

## Chest Drainage Program

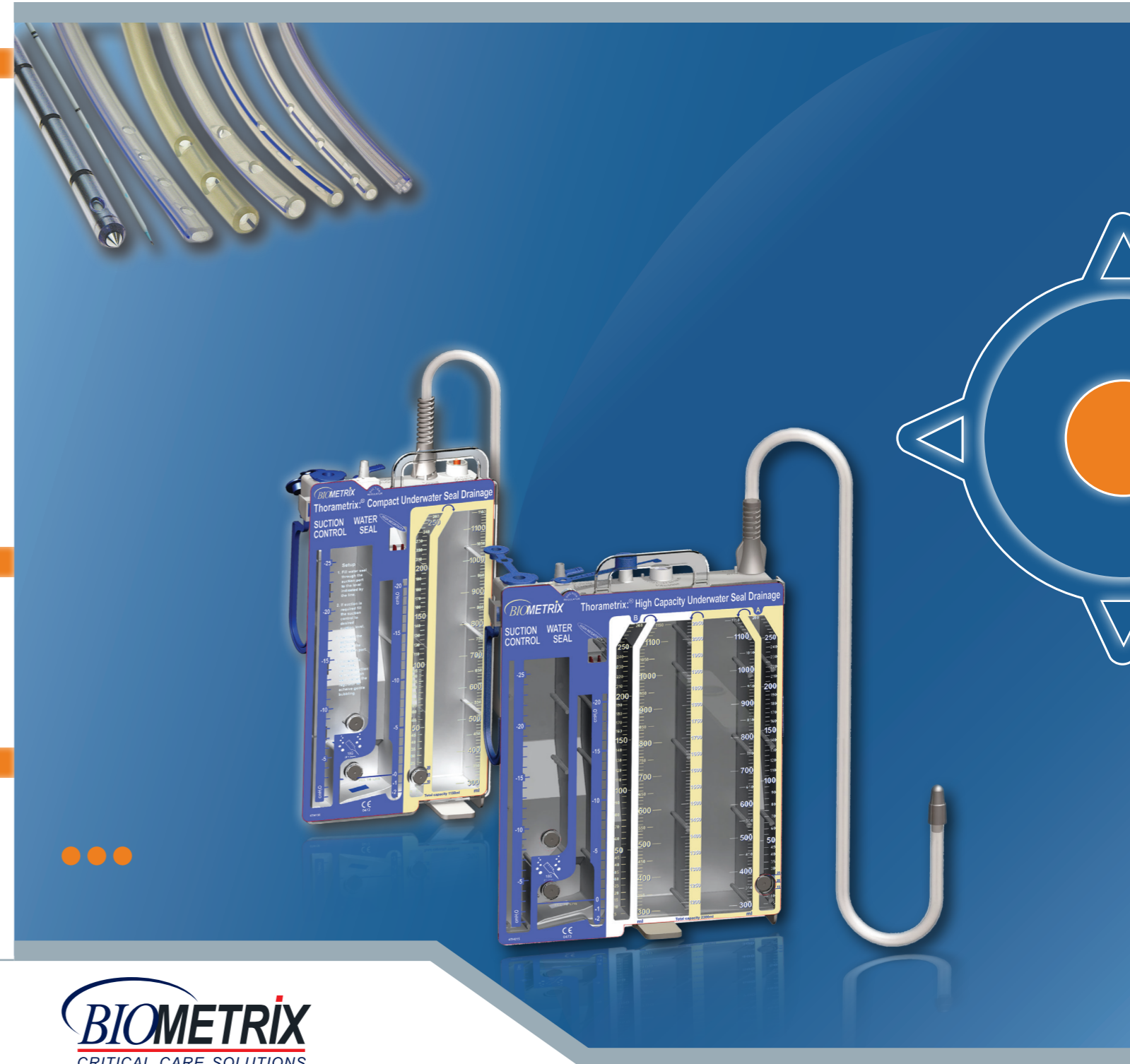

# ... THORAMETRIX™ - Chest Drainage Systems

The Thorametrix™ system is a post operative chest drainage system that is used for draining air or liquids from the pleural cavity. It presents a simple and convenient drainage system for maintaining the required amount of vacuum in patients' chests while controlling the level of suction.


With its clear chambers for accurate fluid loss assessment, easy-to-use handle bed hangers, floor stands and pre-connected drainage tubing, the Thorametrix offers a good value for physicians and patients.

## The Thorametrix chest drainage system features the following benefits:

- Suction Regulator – convenient and easy-to-use adjustment tab for decreasing the bubbling noise and water evaporation from the suction control chamber
- A special secure cap on port B ensures patient safety in case of single drain connection
- Positive pressure automatic relief valve avoids the occurrence of tension pneumothorax
- Double symmetrical chest drainage enables parallel drainage from two different sites of the pleural cavity (in high capacity systems)
- Anti-reflux valve for indicating abnormal high negativity and preventing water loss
- High negativity manual relief valve counteracts pressures generated by patient and prevents loss of seal due to siphoning into collection chamber



**Compact**  
*Ideal for newborn, pediatric, trauma and spontaneous Pneumothorax.*  
*Total capacity: 1150ml*

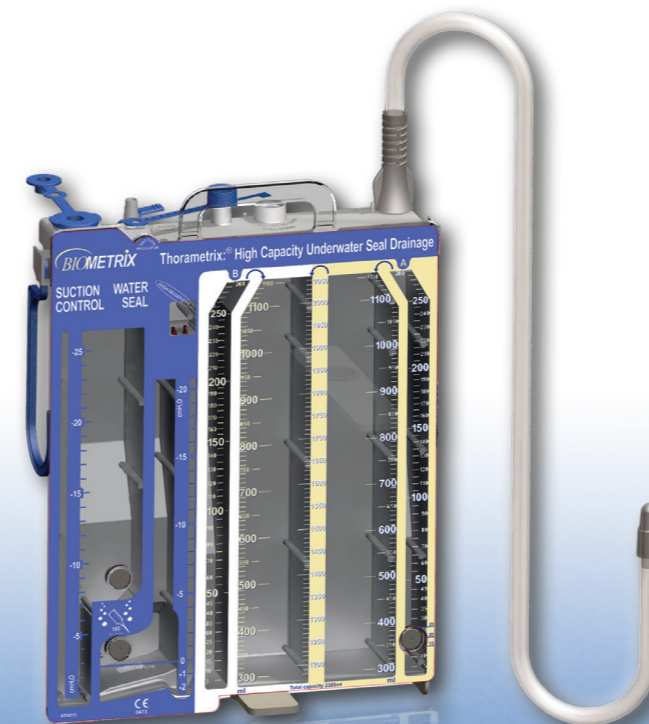

**High Capacity**  
*Ideal for intensive care and cardio-thoracic units.*  
*Total capacity: 2300ml*